# PremEco<sup>®</sup> Line

## Gießküvette

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- |                               |                                  |
|-------------------------------|----------------------------------|
| (en) Instructions for use     | (lt) Naudojimo instrukcija       |
| (fr) Notice d'utilisation     | (lv) Lietošanas informācija      |
| (it) Istruzioni per l'uso     | (et) Kasutusjuhend               |
| (es) Instrucciones de uso     | (hu) Használati utasítás         |
| (pt) Instruções de utilização | (cs) Návod k použití             |
| (pl) Instrukcja użycia        | (sl) Navodila za uporabo         |
| (nl) Gebruiksaanwijzing       | (sk) Návod na používanie         |
| (sv) Bruksanvisning           | (el) Οδηγίες χρήσης              |
| (da) Brugsanvisningen         | (tr) Kullanım kılavuzu           |
| (fi) Käyttöohje               | (uk) Інструкція для застосування |

# PremEco® Line Gießkuvette

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# PremEco® Line Casting Flask

*A robust, non-distorting precision duplicating flask for the resin casting technique for fabricating partial and total resin and model cast prostheses, for double-crown and cover-denture prostheses and for bite raisers of all kinds.*

## **(en)** Instructions for use, please read carefully!

### **Areas of Application**

For the rational completion in the field of polymer prosthetics based on gel and hydrocolloid embedding requiring time-saving preparation and finishing after polymerisation.

Suitable for use with methylmethacrylate-based prosthetic basic resins.

### **Handling**

- Opening and closing the flask using the three closure screws.
- Observe the markings on the side when assembling.
- Agar-agar-based hydrocolloid is recommended for filling the flask.


The PremEco® Line Casting Flask is a module from the Karl-Heinz Körholz PremEco® Line fabrication system. For further information please visit our Internet site at [www.merz-dental.de](http://www.merz-dental.de).

### **Supply Units:**

- 1 casting flask, normal (Ø 110 mm) incl. 1 screwdriver and 1 channel cutter
- 1 casting flask, large (Ø 125 mm) incl. 1 screwdriver and 1 channel cutter
- 1 channel cutter

### **Art.-No.:**

- 1030078
- 1030079
- 1030080
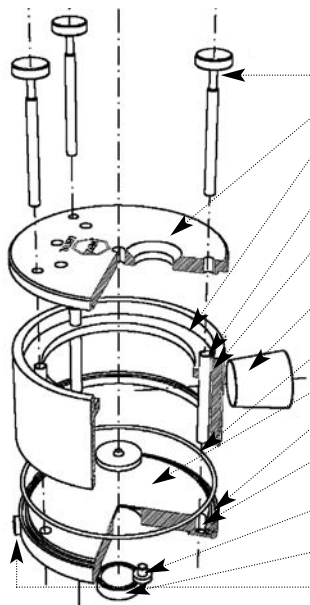
### **The flask is constructed as shown below:**

Three knurled-head screws

Robust flask cover

Stabilizing bar for duplicating gel when opening the flask  
occlusion protection against polymerisation pressure

Fitted sleeves for exact positioning of cover,  
upper part and bottom

Markings on the housing and flask bottom  
facilitate assembly

2 large sealing plugs close the fill openings for resin and provide  
an adequate shrinkage reservoir

O-ring ensures seal of the flask and thus polymerisation  
pressure on one side on the prosthesis base

Aluminium flask floor for better temperature conductivity

Groove for O-ring

Fitted bore holes for establishing anti-twist safeguard in the  
model flask bottom

3 stainless steel adjusting feet for controlled, uniform cooling  
of the gel

Removable magnet for securing the model

Stainless steel feet ensure a secure position when filling

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**Merz Dental GmbH**