Baltic Denture System

Digital solution for dentures

BDKEY[®] Plane & BDKEY[®] Fin

(de) Gebrauchsinformation

(en) Instructions for use

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(sv) Bruksanvisning

(da) Brugsanvisning

(no) Bruksanvisning

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(lt) Informacinis lapelis

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(et) Kasutusteave

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(cs) Informace k použití

(sl) Navodila za uporabo

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BDKEY[®] Plane & BDKEY[®] Fin

en Instructions for use, please read carefully!

Technical applicability


The BDKEY[®] Plane & BDKEY[®] Fin set serves as an optional tool for visualization of the occlusal plane and centre of face. It is intended for use together with the self-printable BDKEY[®] Innerplane.

BDKEY[®] Plane & BDKEY[®] Fin consists of

- 1 x BDKEY[®] Plane
- 1 x BDKEY[®] Fin

Intended users

Dental professionals

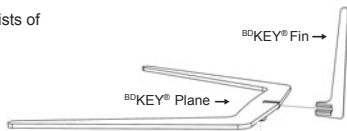Application instructions

- The BDKEY[®] Fin can be attached at the designated spot, centred on the BDKEY[®] Plane.
- The BDKEY[®] Plane, BDKEY[®] Fin & BDKEY[®] Innerplane can only be attached in one direction.
- The BDKEY[®] Innerplane may be optionally stored as a dataset at www.merz-dental.de.
- The BDKEY[®] in sizes S, M and L may be optionally stored as a dataset at www.merz-dental.de.

Storage

Protect against direct sunlight..

The BDKEY[®] Plane & BDKEY[®] Fin are cleaned and visually inspected before use as follows.

Cleaning and visual inspection

- Rinse off soiling under running water; take special care given the presence of cavities.
- Cleaning can be carried out with a soft brush.
- No additional cleaning agents should be used.
- After cleaning, the products should be visually inspected for residues or damage.

For detailed information please see the applicable Material Safety Data Sheet, available upon request or as download from www.merz-dental.de

Disposal

The BDKEY[®] Plane & BDKEY[®] Fin components are water-insoluble, inactive, do not pose a threat to groundwater and may therefore be disposed of with normal household waste.

The product properties are based on compliance with the instructions for use

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Merz Dental GmbH